

The Kinetics of a Diffusion-Controlled Reaction

Dustin Wheeler


March 28, 2025

The objective of this lab is to determine the rate constant and collision diameter of a diffusion-controlled reaction (photoexcited anthracene with carbon tetrabromide) using the technique of fluorescence quenching.¹

¹ Transcribed (with corrections) from Halpern [1].

Introduction

In this experiment, you will study a very fast bimolecular reaction between two different species in solution. Let us assume that the reactants, A and B, which are electrically neutral, undergo independent, random motion in solution. There is a certain probability that A and B will encounter each other at some close distance, R , where R is approximately equal to the sum of their molecular radii. This arrangement is called an *encounter complex*. Because there is a tendency for the solutes A and B to maintain constant random motion, it is inevitable that, once loosely held in this complex, they will subsequently separate from each other unless a chemical reaction (or other definitive process) first links them together or causes them to react. Since this random motion takes place in a “bath” of solvent molecules (assumed to be unreactive with respect to A and B), the separation of the A–B encounter complex will be impeded by the neighboring solvent molecules. This artificial holding-together of the two molecules is called the *cage effect*. The important point is that there is a kinetic competition between a net (thus measurable) reaction between A and B via the collision complex, and the release of A and B from the solvent cage back into the solvent medium where no reaction occurs. We can represent these processes by the following kinetic scheme:

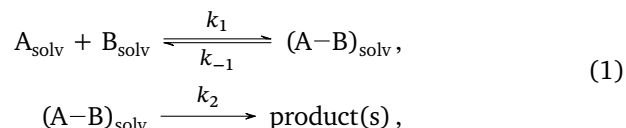

where A_{solv} and B_{solv} represent the solvated A and B molecules and $(\text{A-B})_{\text{solv}}$ is the solvated encounter complex. The rate constants for these elementary steps are denoted as k_1 , for the bimolecular formation of the encounter complex (or diffusion into the solvent cage); k_{-1} , for the unimolecular dissociation of the complex (or diffusion out of the solvent cage); and k_2 for the “unimolecular” reaction between A and B in the complex to form the product(s).

Because we assume that the encounter complex undergoes rapid deactivation, either by dissociation or via reaction, we can employ the

steady-state approximation, according to which the net formation rate of $(A-B)_{\text{solv}}$ is zero. Thus

$$\frac{d[A-B]}{dt} = k_1[A][B] - (k_{-1} + k_2)[A-B] = 0,$$

and we can approximate the steady-state concentration of the encounter complex as

$$[A-B] \simeq \frac{k_1[A][B]}{k_{-1} + k_2}.$$

Now, if we further assume that the reaction rate is given by the slow step, i.e., $k_2[A-B]$, we can express that rate as

$$\text{Rate} \simeq k_2[A-B] = \frac{k_1 k_2}{k_{-1} + k_2} [A][B] = k_{\text{obs}} [A][B] \quad 1/\text{s}.$$

For convenience, we define the second-order rate coefficient, k_{obs} , as

$$k_{\text{obs}} = \frac{k_1 k_2}{k_{-1} + k_2}. \quad (2)$$

This analysis leads to two limiting cases with respect to k_{obs} :

1. the reaction between A and B is very slow compared with their departure (and separation) from the solvent cage, i.e., $k_{-1} \gg k_2$; and
2. the A–B reaction is much faster than their separation, or $k_2 \gg k_{-1}$.

The first case describes a reaction that is under *chemical* control, with $k_{\text{obs}} \simeq k_1 k_2 / k_{-1}$, and the second pertains to a *diffusion-controlled* reaction for which $k_{\text{obs}} \simeq k_1$. The latter situation is considered in this experiment. More details of the reaction, which involves fluorescence quenching, will be described after a discussion of the relevant theoretical ideas.²

Diffusional Mass Transport

The basic issue in a diffusion-controlled reaction concerns the dynamics of mass transport in the condensed phase. The fundamental equations describing mass transport are embodied by Fick's laws of diffusion. Fick's first law says that the number of molecules, dn , diffusing through a unit area, A , per unit time, dt , in the x direction (i.e., the flux, J_x) is proportional to the concentration gradient at that point, i.e.,

$$J_x \equiv \frac{1}{A} \frac{dn}{dt} = -D \frac{dC}{dx}, \quad (3)$$

where the proportionality constant D is called the *diffusion coefficient*. The units of D are cm^2/s . The minus sign in eq. (3) indicates that transport goes against the concentration gradient (i.e., from high to low concentration values).

² Note from the discussion that, when the actual reaction occurs fast enough within the solvent cage (i.e. when k_2 is very large compared to k_{-1}), $k_{\text{obs}} \approx k_1$. From the definition contained in eq. (1), k_1 is the rate constant for the diffusion of the reactants through the solvent. Therefore, this condition means that this diffusion step is actually rate limiting for the overall reaction. Once the reactants are close enough to one another, they will react almost instantaneously. This also implies that such a reaction would have essentially no activation barrier.

Fick's second law states that the change in the concentration of molecules, dC , diffusing across an infinitesimally thin plane per unit time, dt , is proportional to the gradient of the flux:

$$\left(\frac{dC}{dt}\right)_x = -\left(\frac{dJ_x}{dx}\right)_t. \quad (4)$$

Again, the minus sign ensures that the concentration increases in time in response to a flux that decreases with increasing x . If we assume D to be independent of x , we may substitute the expression for J_x from Fick's first law into the second to give

$$\left(\frac{dC}{dt}\right)_x = D\left(\frac{d^2C}{dx^2}\right)_t. \quad (5)$$

Because we are concerned with a three-dimensional and isotropic space, i.e., equal forces in all directions, we can write Fick's second law in a more general form as

$$\left(\frac{dC}{dt}\right)_{x,y,z} = D(\nabla^2 C)_t, \quad (6)$$

where ∇^2 is the Laplacian operator $\partial^2/\partial x^2 + \partial^2/\partial y^2 + \partial^2/\partial z^2$.

The solution of Fick's second law is an equation that expresses concentration as a function of *space* (i.e., distance), r , and *time*, t . We are interested in solving eq. (6) with the boundary conditions

$$\begin{aligned} C(r, 0) &= C_0, \\ C(\infty, t) &= C_0, \\ C(r = R, t) &= 0, \end{aligned}$$

where $r = \sqrt{x^2 + y^2 + z^2}$, and C_0 denotes the *bulk* concentration.

The first boundary condition states that, initially, the bulk concentration prevails throughout the system; the second says that, even after the reaction starts ($t > 0$), the concentration very far away from the reactant, i.e., the bulk concentration, C_0 , is constant; and the third indicates that the reactant concentration is zero at a distance equal to the sum of the reactant collision radii. Note that the diffusion coefficient, D , contained in eqs. (3), (5) and (6) is the sum of the individual diffusion coefficients of the reactants, $D_A + D_B$. This expression essentially allows the motion of one of the reactants to be considered relative to the other.

The solution of Fick's second law for this set of boundary conditions (first rendered by Smoluchowski in 1917) yields, for the space-time dependence of C ,

$$C(r, t) = C_0 \left[1 - \frac{R}{r} \operatorname{erfc} \left(\frac{r - R}{2\sqrt{Dt}} \right) \right], \quad (7)$$

in which $\operatorname{erfc}(x)$ is the *complementary error function*,³ $1 - \operatorname{erf}(x)$, defined as

$$\operatorname{erfc}(x) \equiv 1 - \frac{2}{\sqrt{\pi}} \int_0^x \exp(-y^2) dy.$$

The time dependence of the reactant at the reaction boundary $r = R$ becomes expressible in terms of the flux of reactant, J_R , or its rate of transport across a hypothetical spherical surface with radius R . Thus

$$J_R = \frac{4\pi R D N_A C_0}{1000} \left(1 + \frac{R}{\sqrt{\pi D t}} \right) \quad \text{molecules/s}, \quad (8)$$

where N_A is Avogadro's number. The units shown (and factor of 1000) pertain if C_0 is expressed in molarity (mol/L) and R and D are expressed in cm and cm^2/s , respectively.⁴ Equation (8) indicates that the flux (molecules per unit area per second) is actually time-dependent. This time dependence comes from $R/\sqrt{\pi D t}$, which is sometimes referred to as the *transient term*. Physically, the transient term accounts for the fact that initially nearby reactant molecules do not have to diffuse through the bulk medium in order to react. After a short time, however, these nearby reactant molecules are depleted, and the flux approaches a constant, or steady-state, value, $4\pi R D N_A C_0/1000$. This can be seen mathematically in that time is represented in eq. (8) as \sqrt{t} .

Because J_R represents the rate of passage of one of the reactants (having a bulk concentration C_0 in M) through a spherical reaction surface with the other reactant at the center, the rate of the reaction represented by eq. (1) is merely J_R itself. Hence, the bimolecular, diffusion-controlled rate constant k_1 is J_R/C_0 ,

$$k_1 = \frac{4\pi R D N_A}{1000} \left(1 + \frac{R}{\sqrt{\pi D t}} \right) \quad 1/(\text{M} \cdot \text{s}). \quad (9)$$

This rate constant is not a true *constant*, however, because of the transient term. However, as t becomes large enough,

$$\lim_{t \rightarrow \infty} k_1 = \frac{4\pi R D N_A}{1000} \quad 1/(\text{M} \cdot \text{s}), \quad (10)$$

and eq. (10) can be used to determine the diffusion-controlled rate constant if the mutual diffusion constant ($D_A + D_B$) and collision radii ($R_A + R_B$) are known.⁵ We emphasize again that the units of k_1 presented above are $1/(\text{M} \cdot \text{s})$, if R and D are expressed in cm and cm^2/s , respectively.

Often, however, R and D are not known, and an indirect method is used to estimate k_1 . This approach, developed by Einstein using Stokes's law, allows the diffusion coefficient to be expressed in terms of a bulk property, the solvent viscosity, η . Thus

$$k \simeq \frac{8\bar{R}T}{3\eta} \quad 1/(\text{M} \cdot \text{s}), \quad (11)$$

³ This is available in the SciPy special functions module as `erfc`. Browse the help page for details on the function and its use.


⁴ To convert from molecules to M, we need to account for the number of particles in a mole (Avogadro's number) and convert from our molecular length units (cm) to proper volume units. It's useful to know that $1 \text{ L} = 1 \text{ dm}^3$ and hence molarity can be expressed as mol/dm^3 .

⁵ An important consideration in eqs. (9) and (10) is the time dependence of k_1 . We do not typically think of rate constants as changing in time, but eq. (9) contains a *transient term*, $R/\sqrt{\pi D t}$. Remember that the diffusion, or movement of the molecules through the solvent, is the slow step of this reaction and is rate determining. However, at the very beginning there will be some fraction of molecules that are already close enough together to react without needing to diffuse. Therefore, the apparent rate of the reaction will be faster initially. As these pairs react, this population is depleted, leaving behind reactants that are farther apart and thus need to undergo diffusion. Thus, the transient term dies out and the reaction finally reaches a purely diffusion controlled regime, described by the rate constant in eq. (10).

where $\bar{R} = k_B N_A$ is the ideal gas constant.⁶ Equation (11), which is often referred to as the Stokes-Einstein-Smoluchowski (SES) equation, holds when the reactants A and B are different and their sizes are larger than that of the solvent molecules. The indicated units are obtained if \bar{R} and η are expressed in SI units, such that $\bar{R} = 8.314 \times 10^3 \text{ L} \cdot \text{Pa}/(\text{mol} \cdot \text{K})$ and η is in $\text{Pa} \cdot \text{s}$.

Fluorescence Quenching

The reaction studied in this experiment takes place between an electronically excited molecule and a quencher, a species that removes the electronic excitation. The advantages of using an electronically excited molecule as one of the reactants are: 1) the reaction takes place only when the system is exposed to light; 2) the concentration of electronically excited molecules can be readily followed fluorometrically (via fluorescence detection); and 3) the reaction itself is intrinsically fast—it has no activation energy *per se*, and is limited by the encounter between the excited molecule and the quencher. Many experiments utilize luminescence quenching as a “kinetic” technique. The overall reaction describing fluorescence quenching can be represented as

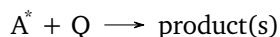


where A^* represents the electronically excited (fluorescent) state of a molecule (e.g., anthracene in this experiment), Q denotes the quencher (e.g., CBr_4) that “extinguishes” the A^* fluorescence, and “product(s)” indicates the species into which the fluorescent molecule is eventually transformed. This de-excitation process may occur through a charge-transfer intermediate ($A^+ - Q^-$) in which electronic charge is transferred to the quencher molecule from the excited-state species. For the system encountered in this experiment, electron transfer from A^* to the quencher occurs very rapidly. In general, it is possible that the consequence of fluorescence quenching is the return of A to its electronic ground state without any change in structure, i.e., $A + Q$.

In the case of CBr_4 , the quenching mechanism involves the electron transfer from the excited state, A^* , to CBr_4 , which is more electronegative than A^* . A reverse electron transfer then restores the two molecules to their neutral species, A being in the ground state.

The bimolecular quenching reaction shown above, however, competes with the intrinsic first-order fluorescence decay of A^* and thus reduces the probability that the A^* species will emit a photon. A more

⁶ Note the distinction between \bar{R} and R . \bar{R} appears *only* in this equation and refers to the gas constant. Everywhere else, R is the summed “collision radii” of the interacting species. It is this collision radius that we are trying to find in this experiment, along with the rate constant k_1 . Also, keep in mind that eq. (11) represents an *alternative* method for estimating the diffusion rate constant, k_1 , based on bulk parameters of the solvent, rather than detailed knowledge of the molecular parameters R and D .

complete scheme is thus

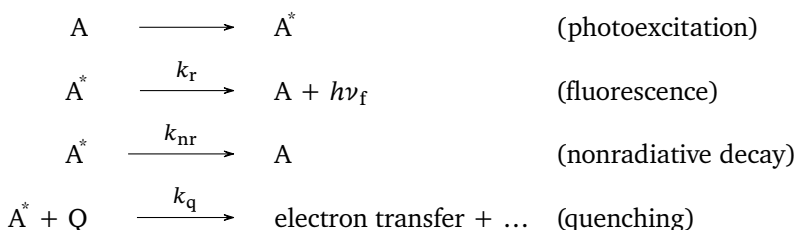

If a solution containing A and Q is irradiated with a *steady-state* light source (at a wavelength where A absorbs), a very small time-independent concentration of A^* is produced. The probability that A^* will fluoresce, P_f , is the ratio of the fluorescence (or radiative) rate constant to the sum of all the rate constants that deplete A^* , namely,

$$P_f = \frac{k_r}{k_r + k_{nr} + k_q[\text{Q}]} \quad (12)$$

where k_r , k_{nr} , and $k_q[\text{Q}]$ are, respectively, the radiative, nonradiative, and (pseudo-first-order) quenching rate constants. Since $[\text{Q}] \gg [\text{A}^*]$, the diffusion-controlled bimolecular quenching step becomes pseudo-first-order. If the amount of light absorbed by A is unaffected by [Q], the fluorescence *intensity*, I , measured by the detector (i.e., the instrument response) is proportional to P_f . In the absence of quencher, $I = I_0 = k_r/(k_r + k_{nr})$, and the ratio of fluorescence intensities of unquenched to quenched samples of A becomes

$$\frac{I_0}{I} = 1 + \frac{k_q[\text{Q}]}{k_r + k_{nr}}. \quad (13)$$

Identifying $1/(k_r + k_{nr})$ as the reciprocal of the *fluorescence lifetime* of A^* in the absence of quencher (τ_0), we can write eq. (13) as

$$\frac{I_0}{I} - 1 = k_q\tau_0[\text{Q}]. \quad (14)$$

Equation (14) is known as the steady-state Stern-Volmer (S-V) relation as applied to fluorescence quenching. It predicts that a plot of $I_0/I - 1$ vs. $[\text{Q}]$ should be linear, having an intercept of zero and a slope equal to $k_q\tau_0$. The latter is called the Stern-Volmer constant, K_{SV} (units of $1/M$). If fluorescence quenching is diffusion-controlled (which is the case in the anthracene/ CBr_4 system), k_q can be identified with k_1 (see eqs. (1), (2) and (9) to (11)).⁷

In the absence of the transient effect, which is a type of “static” quenching in that some A^* molecules are quenched by nearby Q molecules that do not have to undergo transport fully through the solvent medium in order to approach A^* within a distance R , the S-V plot should be linear, and independent knowledge of τ_0 allows us to determine k_q . Alternatively, if the process is known to be diffusion-controlled, and if k_1

⁷ The first Stern-Volmer Equation eq. (14) is generally applicable because it makes no assumptions about the quenching rate constant k_q . Since we are dealing with a diffusion-controlled reaction, we can identify k_q with the diffusion-controlled rate constant, k_1 , from eqs. (9) to (11). It should therefore be straightforward to measure the fluorescence emission of some analyte A as we change the concentration of the quencher Q, and fit the data to eq. (14) to obtain a value for k_1 . However, this still leaves us with a question: Does the k_1 obtained in this manner correspond to eq. (9), eq. (10), or eq. (11)?

Equations (10) and (11) do not depend on the time, t , and therefore correspond to the steady-state regime where the reaction is fully diffusion-controlled. When combined with eq. (14), they describe a linear Stern-Volmer plot in which the slope is equal to $k_1\tau_0$.

Equation (9), however, includes a time dependence of k_1 . The expression cannot simply be inserted into eq. (14) without carrying this time dependence along with it. Since we are not measuring the quenching as a function of time, this doesn't help us. Instead, the Stern-Volmer Equation can be re-derived, taking into account the transient term, to obtain eq. (15).

is determined from the SES (or some other) equation (e.g., eq. (10)), the fluorescence lifetime of A^* can be estimated from the measured K_{SV} value. This is a case in which a dynamic property (a rate constant) can be obtained from a static experiment (I_0/I measurements).

The S-V relation can be expressed in terms of the diffusion-controlled rate constant that contains the transient term, e.g., eq. (9). The time dependence was treated by Ware and Novros [2], who applied eq. (9) to the condition pertinent to this experiment, namely, that the system is considered to be suddenly exposed to a steady-state excitation source, called a step function (see fig. 1 and Further Readings).

In this case, the S-V equation becomes more complicated:

$$\frac{I_0}{I} = \frac{1 + (4\pi RDN_A/1000)\tau_0[Q]}{Y}, \quad (15)$$

where

$$Y = 1 - \frac{b\sqrt{\pi}}{\sqrt{a}} \exp\left[\left(\frac{b^2}{a}\right)\right] \operatorname{erfc}\left(\frac{b}{\sqrt{a}}\right), \quad (16)$$

in which

$$a = \frac{1}{\tau_0} + \frac{4\pi RDN_A}{1000}[Q] \quad 1/s,$$

$$b = \frac{4\sqrt{\pi DR^2 N_A}}{1000}[Q] \quad \sqrt{1/s}.$$

If $Y = 1$, eq. (15) reduces to eq. (14),⁸ in which case $k_q = k_1$ from eq. (10).⁹

The Stern-Volmer plot based on eq. (15) should be curved upward. If the mutual diffusion coefficient of the reacting pair is known, the S-V plot of experimental data can be fit using eq. (15) with R used as a “fitting” parameter. For a satisfactory match between experiment and theory, however, the value of R must be physically reasonable. That is, R should be approximately equal to the sum of the molecular diameters of A (actually A^*) and Q .

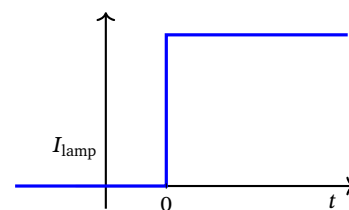


Figure 1: A “step excitation” function (known in mathematics as a “Heaviside function”). At time $t = 0$, absorbing radiation is suddenly “switched on.”

⁸ Equation (15) describes a curved line due to the transient effects that take place very soon after the light source is switched on. However, this may cause some confusion, because our Stern-Volmer plots do not involve a time variable. How can the transient effects account for this curvature, if there is no explicit time dependence? We touched on the answer to this before, but it is worth emphasizing again. The transient term arises because when the light source is first switched on there are some A/Q pairs that are close enough to react without diffusing very far. Since the overall rate of the reaction can be thought of as an average of the rates for each individual molecular collision, the faster reactions between these pairs will skew the average towards faster rates. This statistical effect will be much more significant when there is a very high concentration of one or both components, and almost negligible when the concentrations are very low. Since the slope of the Stern-Volmer plot is proportional to the rate constant, k_1 , this means that the slope will be large (faster rates) at high quencher concentrations, $[Q]$, and small (slower rates) at low $[Q]$. We will explore how these conditions relate to the other Stern-Volmer equation (eq. (14)) in the Data Analysis section below.

⁹ Verify this for yourself!

Safety Precautions

- Always wear safety glasses or goggles; these glasses should block ultraviolet light. Ordinary plastic safety goggles or glasses may not be effective in absorbing all the ultraviolet radiation. Check with your instructor.
- Do not allow solid anthracene to come in contact with the skin. Wear gloves when handling anthracene and CBr_4 .
- Ozone is sometimes produced by ultraviolet light sources. If you detect this gas, which has a sharp, slightly acrid odor, notify your instructor and immediately shut off the source. Increase ventilation, and leave the room immediately.
- The compressed-gas cylinder used in deaerating should be securely strapped to a firm support. The delivery pressure should never exceed a few (e.g., 5) psig.
- The experiment should be performed in an open, well-ventilated laboratory.

Procedure

Required Equipment

- (10) 10 mL volumetric flasks
 - (1) 10 mL graduated pipet
 - (1) fluorescence cell
 - dry nitrogen gas cylinder
1. Make sure the fluorimeter is turned on and warmed up before taking any spectra. Once the instrument is warmed up, you'll want to check the baseline zero level prior to collecting any data.
 2. Using freshly sublimed or purified anthracene and CBr_4 , prepare 250 mL of $\sim 1.00 \times 10^{-4} \text{ M}$ solution of anthracene (AN) in spectro-metric (or fluorometric) quality *n*-hexane. Using this solution as "solvent," prepare 25 mL of $\sim 1.50 \times 10^{-2} \text{ M}$ "stock" solution of CBr_4 .
 3. Prepare at least eight dilutions of the AN/ CBr_4 stock solution using the AN solution as solvent. The most dilute should be 5 %, and the rest should be nearly evenly spaced up to the stock-solution value. A volume of 10 mL of each solution is appropriate. Include 0 % and 100 % samples for uniformity. Label these samples.

4. Starting with the 0% sample (AN only), introduce the solution into a fluorescence cell and deaerate with dry N_2 for about 2 min. Avoid using an excessive flow rate that will cause the solution to splatter from the cell or to otherwise result in undue solvent evaporation. Promptly stopper the cell.
5. Record the full fluorescence spectrum using the instrumental conditions previously outlined by the instructor.¹⁰ Make sure the spectrum record is labeled.
6. Using exactly the same settings as in step 4, measure the fluorescence intensity of the most dilute and subsequently more concentrated AN/CBr₄ solutions. You must deaerate each sample in a consistent procedure (i.e., identical bubbling time). If you have any doubt as to the steadiness of the exciting source, switch back to the pure AN sample and compare its maximum intensity with that of the first sample run.
7. If your spectrometer has a manual control for the shutter, make sure it is *closed* until you're ready to start each measurement, at which point you should open it. When the measurement is complete, close the shutter to prepare the spectrometer for the next sample.
8. Keep in mind the following procedural points:
 - (a) The temperature of the sample should be held constant within 0.5 °C to 1.0 °C;
 - (b) the exposure times of all of the samples (especially those rich in CBr₄) to the excitation source (and even room light) should be minimized; and
 - (c) samples should be stored in the dark until they are used.

Parameter	Value
λ_{exc}	283 nm
λ_{det}	350 nm to 550 nm
Scan speed	300 nm/min
Slit width	2.5 nm

Data Analysis

1. Tabulate your data as I (the anthracene fluorescence at 398 or 399 nm) and $[Q]$ (the CBr₄ concentration, in M).
2. Plot a graph of I_0/I vs. $[Q]$. Note that I_0 is simply the value of the unquenched fluorescence peak (when $[Q] = 0$). This is your experimental quenching curve or Stern-Volmer plot, and we will now see how it corresponds to the two limits of the quenching model (with and without the transient term), as well as the SES equation.
3. Begin by considering the Stokes-Einstein-Smoluchowski (SES) equation, eq. (11). Use the gas constant, \bar{R} , the temperature (assume 25 °C), and the viscosity of *n*-hexane to calculate k_1 . Use the proper

You can then find the wavelength for the maximum fluorescence and grab the intensity at that wavelength for each sample.

units as described in the manual. Once you have k_1 , insert it into eq. (10) to calculate the collision radius R (don't mix up \bar{R} and R).¹¹ What is the collision radius calculated this way (including the unit)?¹²

4. Using the k_1 you just calculated from the SES equation, superimpose the straight line obtained from eq. (14) onto your Stern-Volmer plot. You'll need the lifetime, τ_0 .¹³ Describe the relationship between this line and the experimental data points, if any.
5. Now we will consider the Stern-Volmer plot obtained when we ignore the transient term. Remember that this condition holds at low concentrations. Therefore we will be looking for the tangent to the experimental curve near the y-axis. Take a few of the experimental data points close to $[Q] = 0$, fit a trend line to them, and show the equation. The points you include for this should all be nearly in a straight line. If you take too many, they will begin to show the curvature and your slope will be incorrect. From the slope of this tangent line, calculate an alternate value for k_1 . You'll need the lifetime from the manual again. From that value of k_1 , calculate another value for the collision radius, R , using eq. (10). How do these values of k_1 and R compare to the values calculated from the SES equation in step 3?
6. Now we can examine the entire data set, which includes effects due to the transient-term, in order to extract a third value of R directly. We will do this by fitting our data to eq. (15) using a "nonlinear regression analysis" with R as the only parameter.¹⁴ Just make sure you realize that the parameter R enters into eq. (15) directly in the numerator, as well as through the denominator Y due to the R dependence of a and b . Thus, you need to make sure that all of these parts of the equation refer to the same variable for R so that it is varied consistently throughout. Make sure that you use $\text{erfc}()$ and not $\text{erf}()$. You can start with a guess of 10 \AA for the radius R and have $\text{curve_fit}()$ solve for the best fit of your line. Superimpose the theoretical curve corresponding to this R value onto your data points.

¹¹ To do this, you'll need to calculate the diffusion coefficient, D , for n -hexane from the D known for n -heptane. Dymond and Øye [3] gives values for the viscosity of various n -alkanes. Ware and Novros [2] reported a value of $D = 4.35 \times 10^{-5} \text{ cm}^2/\text{s}$ in n -heptane at 25°C . To obtain D for n -hexane at some other temperature, assume that D is proportional to T/η . They also reported a value of $\tau_0 = 5.52 \times 10^{-9} \text{ s}$. You can assume this to be temperature independent.

¹² **Keep track of your units, especially moles vs. molecules.**

¹³ To do this, you'll need to generate a set of x -values for your function, then apply your function to those values. The help page for $\text{erfc}(x)$ linked above shows how to do this for that particular function.


¹⁴ To perform this step with Python, you will need to use the `curve_fit()` module from the SciPy library. Instructions and an example can be found in the SciPy Documentation under [scipy.optimize.curve_fit](https://docs.scipy.org/doc/scipy/reference/optimize/curve_fit.html).

Questions and Further Thoughts

1. Why is the anthracene fluorescence quenching diffusion-controlled as opposed to chemical-controlled? I.e., justify why the reaction between AN^* and CBr_4 is very fast.
2. What experiment(s) could you perform to determine whether a

stable photoproduct is formed as a result of the fluorescence quenching?

3. What could you do to obtain supporting evidence for the existence of a charge-transfer (or ion-pair) intermediate, e.g.,

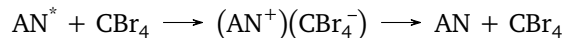

in the quenching process?

4. Verify that if $Y = 1$ in eq. (15), $k_q \rightarrow k_1$ (see eqs. (10) and (14)).
5. One of the undesirable complications associated with this type of experiment is the potential presence of *ground-state* complexes between AN and CBr_4 . If such a complex competes for light absorption with uncomplexed, free anthracene, it can cause a decrease in fluorescence intensity without the existence of diffusion between AN^* and CBr_4 . This is sometimes called *static quenching*. For this reason, it is desirable to work with quencher concentrations as low as possible. What experiment could you do to obtain evidence of a ground-state complex? How could you determine the equilibrium constant of such a complex?
6. What effect would a polar solvent, e.g., acetonitrile (CH_3CN), have on the ground- and excited-state interactions between AN and CBr_4 ?
7. If you wanted to do a fluorescence quenching experiment that demonstrated (and thus maximized) the transient effect in diffusional processes, what qualities of solvent viscosity and fluorescence probe lifetime would you seek (i.e., high or low viscosity; long or short fluorescence lifetime)?
8. The lowering of dissolved oxygen by bubbling the solution with dry N_2 gas (deaerating) is an application of Henry's law. Explain how deaeration works in this context.

Lab Report Guidelines

Your lab report should consist of the following parts:

Title, Author and Date

Experimental Procedure This should be a very brief general outline of the procedure, written out as a paragraph or two. Give the make and model for any major instruments you used, as well as any important settings. For fluorescence spectroscopy, this especially means the excitation wavelength and slit widths.

Results and Discussion This should include the following:

1. A Stern-Volmer graph consisting of:
 - (a) Your experimental data points
 - (b) The straight line from step 4 corresponding to the SES
 - (c) The tangent line from step 5 corresponding to the steady-state regime near $[\text{CBr}_4] = 0$.
 - (d) The theoretical curve from step 6 corresponding to Eq. 15.
 These lines and points should all be styled uniquely to distinguish them, and the equations displayed on or near the chart. Label all axes appropriately (with units).
2. A table of k_1 and R values calculated in the different ways discussed. You don't calculate k_1 from eq. (15), so you will only have two different values for k_1 , but three different values for R . Also give the value of D that you calculated for *n*-heptane, with proper units (this can be given in a footnote to the table, or just the body of the discussion itself).
3. Estimate the errors in each of these tabulated values as best you can, and include them in the table. Do they agree to within the errors or uncertainty? Which value of R should be the most accurate? Which do you trust the most based on your handling of the experiment and analysis?
4. What are the approximate radii of the anthracene and CBr_4 molecules? Since the observed collision radius, R , should be the sum of these two radii, do your experimental values make sense?
5. Answer any four of the questions in the section "Questions and Further Thoughts". This should not be a separate section, but should instead be included organically in the discussion as a way of filling it out.¹⁵

References

Appendix At the very end of your report, include examples of any calculations that you did by hand. Provide digital copies of the Excel (or other) files that you used to generate your graphs.

You do *not* need to include uncertainty calculations for this lab.

References

- [1] Arthur M Halpern. *Experimental Physical Chemistry: A Laboratory Textbook*. Prentice Hall, New Jersey, 1997.

¹⁵ Do **not** submit a lab report with a section labeled "Answers to Further Questions", or some such nonsense... Your goal is to present a lab report akin to those presented in scientific journals. If you're not familiar with the style, skim the latest issue of *J. Chem. Phys.*

- [2] William R. Ware and Joel S. Novros. “Kinetics of Diffusion-Controlled Reactions. An Experimental Test of the Theory as Applied to Fluorescence Quenching”. In: *The Journal of Physical Chemistry* 70.10 (1966), pp. 3246–3253. DOI: [10.1021/j100882a038](https://doi.org/10.1021/j100882a038).
- [3] J. H. Dymond and H. A. Øye. “Viscosity of Selected Liquid *n*-Alkanes”. In: *Journal of Physical and Chemical Reference Data* 23.1 (1994), pp. 41–53. DOI: [10.1063/1.555943](https://doi.org/10.1063/1.555943).

Further Reading

- Joseph R Lakowicz. *Principles of Fluorescence Spectroscopy*. Springer Science+Business Media, 2006.
- *Modern Molecular Photochemistry*. 1st ed. Benjamin-Cummings, Menlo Park, CA, 1978, p. 311.
- RM Noyes. “Effects of diffusion rates on chemical kinetics”. In: *Progress in Reaction Kinetics*. Vol. 1. Pergamon Press, 1961, pp. 129–160.
- JB Birk, ed. *Photophysics of Aromatic Molecules*. John Wiley Sons Inc., New York, 1970, p. 518.
- JH Noggle. *Physical Chemistry*. 3rd ed. Harper Collins, New York, 1996, p. 480.
- PW Atkins. *Physical Chemistry*. 5th ed. Oxford University Press, Oxford, 1994, p. 795.
- RA Alberty and RJ Silbey. *Physical Chemistry*. 1st ed. John Wiley Sons Inc., New York, 1992, pp. 692–696.
- JB Birk. *Organic Molecular Photophysics*. Vol. 1. 1973, p. 403.
- IN Levine. *Physical Chemistry*. 4th ed. McGraw-Hill, New York, 1995, p. 469.